




1285 Corporate Center Drive #150
Eagan, MN 55121 U.S.A.
1-800-624-5662 TOLL FREE IN U.S.A.
1-612-452-1977
1-612-452-1787 FAX

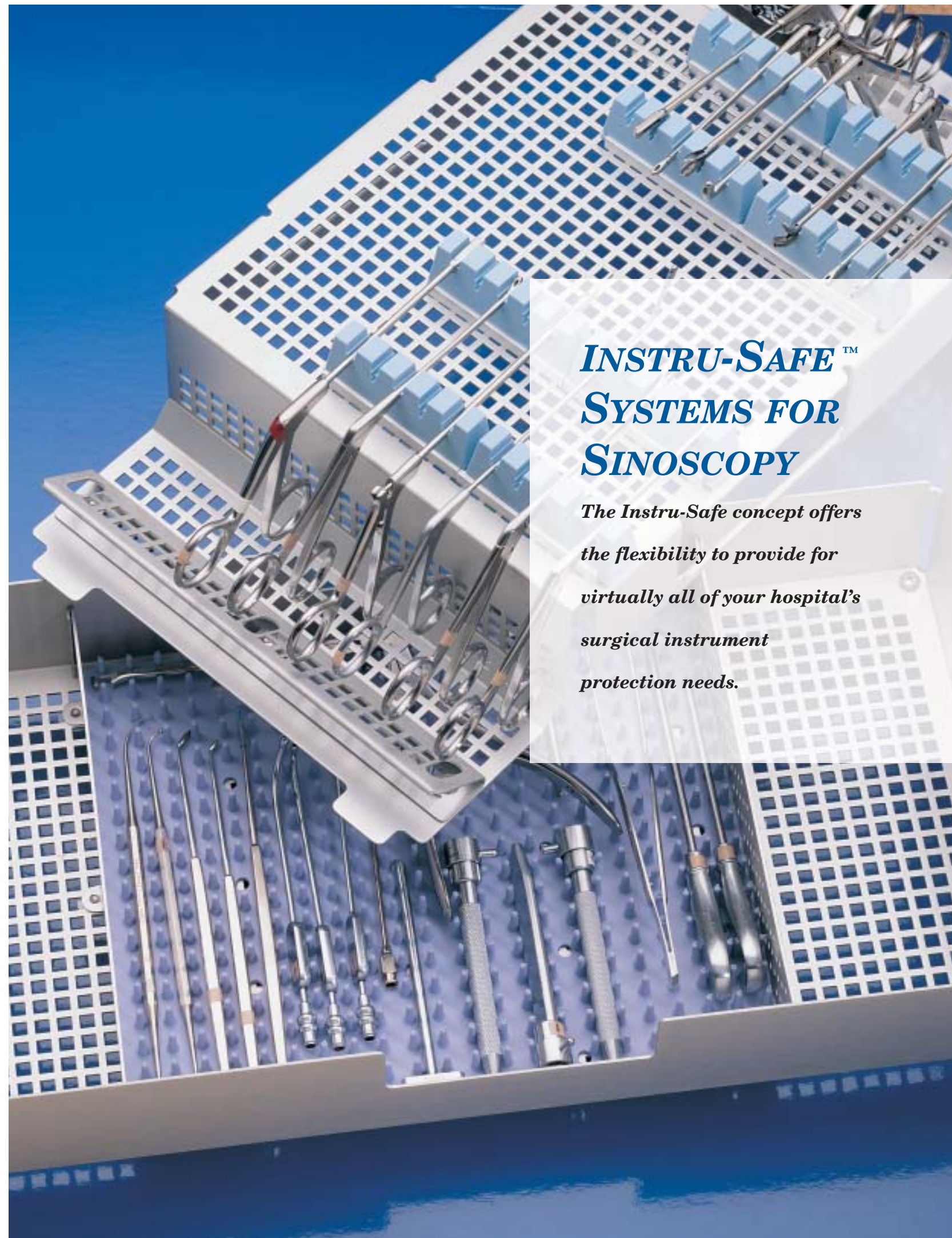

INSTRU-SAFE™ SYSTEMS FOR SINOSCOPY

*The Instru-Safe concept offers
the flexibility to provide for
virtually all of your hospital's
surgical instrument
protection needs.*

INSTRU-SAFE PROTECTIVE TRAYS ARE SUPERIOR BY DESIGN.

Superior Protection

- Reduces instrument repair and replacement costs.
- Keeps instruments in best operating condition.

Superior Organization

- Instruments can be kept together in sets during processing, making instrument accountability easier.

Superior Processing

- Designed to allow sterilant flow and to prevent moisture from accumulating.

Superior Versatility

- Models to accommodate all instruments.
- Usable in any sterile container or wrap.
- Can be sterilized in steam, gas, plasma or disinfectant soak.
- Customized for any application.

Superior Durability

- Anodized aluminum and molded silicone withstand years of use without wearing out or breaking down.


Superior Value

- Long life; pays for itself in reduced repair and replacement.



HOLDS AN ENTIRE SET of instruments normally required for sinoscopy procedures. IN-8810 accommodates up to 20 ring handled instruments plus other associated instruments. Capacity can be increased up to 30 instruments with additional Hold-Its (See accessory line). Can be sterilized in steam, gas, plasma or disinfectant soak. Fits into most full size sterilization container systems or may be wrapped. Stackable.

Order No.	Description	Dimensions
IN-8810(*)	Sinoscopy System	19 3/8" x 9 3/8" x 5 1/2"

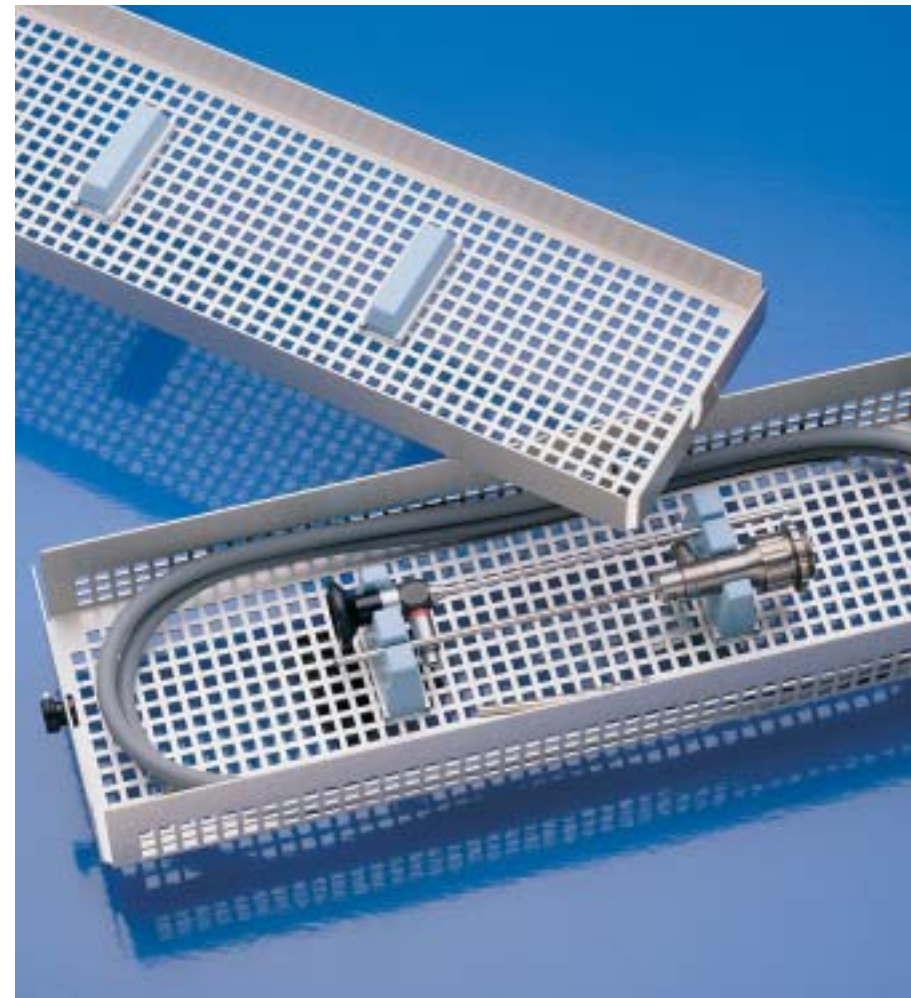

HOLDS ONE OR TWO endoscopes with associated light cables. Accommodates various makes and models of scopes. Can be sterilized in gas or disinfectant soak. Nonstacking.


Order No.	Description	Dimensions
IN-8642	Scope System, Narrow	17 3/4" x 5 3/8" x 2 3/8"

OTHER STANDARD SYSTEMS

MICROMEDICS OFFERS a full line of standard Instru-Safe products to protect and organize surgical instruments.

CUSTOM SYSTEMS

Micromedics can build custom systems to your specific needs. Simply contact your local Micromedics representative, or call Micromedics at 800-624-5662.

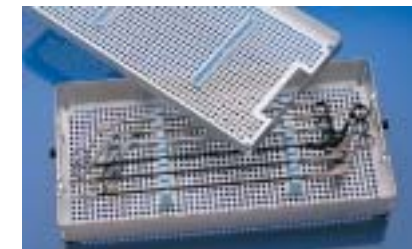

(*) ORDERING INFORMATION

Type of foot provided depends on how system is to be used. Please specify one of the following at time of order:

- AE — For systems that will be used in an Aesculap container
- GN — For systems that will be used in a Genesis container
- MW — For systems that will be used in a Monarch container; or stacked and wrapped
- CF — For systems that will be wrapped only (can not be stacked)



MICROMEDICS OFFERS A VARIETY of accessories which are used to customize systems to your specific requirements. Simply contact your local Micromedics representative, or call Micromedics at 800-624-5662.